



Exploration of xenobiotic metabolism within cell lines used for Tox21 chemical screening

Wei Qu^a, David M. Crizer^a, Michael J. DeVito^a, Suramya Waidyanatha^a, Menghang Xia^b, Keith Houck^c, Stephen S. Ferguson^{a,*}

^a Division of the National Toxicology Program, National Institute of Environmental Health Sciences, Research Triangle Park, NC, USA

^b Division for Pre-Clinical Innovation, National Center for Advancing Translational Sciences, National Institutes of Health, Bethesda, MD, USA

^c National Center for Computational Toxicology, US EPA, Research Triangle Park, NC, USA

ARTICLE INFO

Keywords:

Tox21
Xenobiotic metabolism
Metabolic clearance
Metabolic competence
CYP450
Primary human hepatocytes
Bisphenol

ABSTRACT

The Tox21 Program has investigated thousands of chemicals with high-throughput screening assays using cell-based assays to link thousands of chemicals to individual molecular targets/pathways. However, these systems have been widely criticized for their suspected lack of 'metabolic competence' to bioactivate or detoxify chemical exposures. In this study, 9 cell line backgrounds used in Tox21 assays (i.e., HepG2, HEK293, Hela, HCT116, ME180, CHO-K1, GH3.TRE-Luc, C3H10T1/2 and MCF7) were evaluated via metabolite formation rates, along with metabolic clearance and metabolite profiling for HepG2, HEK293, and MCF-7aroERE, in comparison to pooled donor (50) suspensions of primary human hepatocytes (PHHs). Using prototype clinical drug substrates for CYP1A2, CYP2B6, and CYP3A4/5, extremely low-to-undetectable CYP450 metabolism was observed (24 h), and consistent with their purported 'lack' of metabolic competence. However, for Phase II metabolizing enzymes and metabolic clearance, surprisingly proficient metabolism was observed for bisphenol AF, bisphenol S, and 7-hydroxycoumarin. Here, comparatively low glucuronidation relative to sulfation was observed in contrast to equivalent levels in PHHs. Overall, while a lack of CYP450 metabolism was confirmed in this benchmarking effort, Tox21 cell lines were not 'incompetent' for xenobiotic metabolism, and displayed surprisingly high proficiency for sulfation that rivaled PHHs. These findings have implications for the interpretation of Tox21 assay data, and establish a framework for evaluating of 'metabolic competence' with in vitro models.

1. Introduction

The Toxicology in the 21st Century Program (Tox21) is a federal collaboration among the National Institute of Environmental Health Sciences, National Center for Advancing Translational Sciences, the U.S. Environmental Protection Agency, and the Food and Drug Administration (Collins et al., 2008). Tox21 researchers utilize alternative high-throughput screening methods to evaluate chemical-induced bioactivities (Kavlock et al., 2009; Tice et al., 2013). These approaches have been effective in characterizing thousands of chemicals with respect to

individual molecular targets and pathways, but have often been criticized for their presumed lack of tissue-like differentiation and functionality that includes xenobiotic metabolism 'competence'. Given the central role of liver in human drug metabolism, drug-drug interactions and toxicity, hepatic metabolism is often the focus of these critiques. While hepatic metabolism (i.e., metabolic clearance, metabolite formation rates/profiles, reactive metabolites) is not directly operative in extrahepatic tissues (e.g., endocrine-related), the notion of conferring 'metabolism' to Tox21 assay systems remains a primary focus of researchers working towards new approach methodologies (NAMs) in the

Abbreviations: ABT, aminobenzotriazole; ACTT, American Type Culture Collection; AOPI, acridine orange propidium iodide; BPA, bisphenol A; BPAF, bisphenol AF; BPAF G, bisphenol AF glucuronide; BPAF S, bisphenol AF sulfate; BPS, bisphenol S; BPS G, bisphenol S glucuronide; BPS S, bisphenol S sulfate; BPS 2S, bisphenol S 2 sulfate conjugates; CYPs, cytochromes P450; DMSO, dimethyl sulfoxide; FBS, fetal bovine serum; 7HC, 7-hydroxycoumarin; 7HCG, 7-hydroxycoumarin glucuronide; 7HCS, 7-hydroxycoumarin sulfate; LC-MS, liquid chromatography coupled to mass spectrometry; MCF7, MCF-7aroERE cells; MI, midazolam; NAMs, new approach methods; NTP, National Toxicology Program; OMP, omeprazole; PHHs, primary human hepatocytes; RIF, rifampicin; SD, standard deviation; SULT, sulfotransferase; Tox21, the Toxicology in the 21st Century Program; UGT, UDP-glucuronosyltransferase.

* Corresponding author at: NTP, NIEHS, PO Box 12233, Mail Code K2-17, Research Triangle Park, NC 27709, USA.

E-mail address: stephen.ferguson@nih.gov (S.S. Ferguson).

<https://doi.org/10.1016/j.tiv.2021.105109>

Received 29 October 2020; Received in revised form 2 February 2021; Accepted 5 February 2021

Available online 18 February 2021

0887-2333/© 2021 Published by Elsevier Ltd.

toxicology field to confer bioactivation of parent chemicals to reactive or bioactive metabolites. However, to date ‘metabolic competence’ in Tox21 assay systems has not been directly evaluated in context with established benchmarks of hepatic metabolic competence. Previous studies of metabolism with immortalized cell lines have been limited, but a recent study evaluating of benzo(a)pyrene metabolism in ARE-bla, GR-bla, and AREc32 cells (i.e., engineered reporter cells in HepG2 and HEK293) reported fluorescent-based CYP1A1 enzymatic activities, which reflect fetal liver metabolism (Fischer et al., 2020).

The benchmark for modeling human liver response to xenobiotic exposures in vitro has centered around isolation and culture of primary liver cells (e.g., those isolated from living tissue). Suspensions of primary human hepatocytes (PHHs) are widely recognized as the ‘gold standard’ to quantitatively model human hepatic drug metabolism, as they retain in vivo-expressed enzymes (e.g., cytochromes P450, CYPs) in the initial minutes/hours post isolation from human liver tissue (Hewitt et al., 2007; Guengerich, 1997). Given the myriad of metabolizing enzymes and transporters involved in drug biotransformation and clearance, sentinel drug metabolizing enzyme activities such as CYP3A4/5, CYP1A2 and CYP2B6 are often evaluated to gauge the fidelity, differentiation status, and metabolic proficiency of in vitro liver models including isolated PHHs from donor tissues (Kafert-Kasting et al., 2006). These 3 major CYPs are involved in clearance of >60% of marketed drugs (Gerin et al., 2013; Lee and Kim, 2013). PHHs are routinely used for drug development and environmental toxicology applications to quantitatively estimate metabolic clearance rates and metabolite formation profiles in humans (Wetmore et al., 2012). In addition, confluent cultures of PHHs faithfully reflecting the in vivo-like cobblestone networks of cell-cell interactions observed in vivo have been widely used to model a broad complement of chemical-biological interactions for drug regulatory studies to estimate the potential for drug-drug interactions (e.g., liver enzyme induction, biliary excretion, and time-dependent inhibition) (Kafert-Kasting et al., 2006; Mao et al., 2016).

In addition, CYP450-mediated xenobiotic metabolism, Phase II conjugation via sulfate, glucuronide, glutathione and other facilitated quenching/excretion products play a key role in xenobiotic exposure dynamics. To evaluate Phase II glucuronidation and sulfation proficiency, reference substrates 7-hydroxycoumarin, bisphenol AF (BPAF), and bisphenol S (BPS) were selected for this study. BPAF and BPS represent 2 environmentally-relevant bisphenol A (BPA) alternatives with existing Phase II metabolism and Tox21 assay data (Yamasaki et al., 2003) (Liao et al., 2012) (Song et al., 2012; Wang et al., 2015). Interpretation of their relative potencies is largely presumed a direct reflection of parent chemical exposures based on nominal parent concentrations in the absence of xenobiotic metabolism. An improved understanding the dynamics of in vitro exposures over time that includes Phase II metabolism of Tox21 assays is essential to contextualize and translate bioactivity assay data to humans towards the emergence of next generation risk assessment frameworks.

In this study, we set out to benchmark the ‘metabolic competence’ of the most prevalent cell lines applied in the Tox21 Program (henceforth abbreviated as “Tox21 cell lines”): HepG2, HEK293, Hela, HCT116, ME180, CHO-K1, GH3.TRE-Luc, C3H10T1/2 and MCF-7aroERE (MCF7). Specific activities for major human drug metabolizing enzymes CYP3A4/5, CYP1A2, and CYP2B6 were evaluated using clinically-relevant probe substrates and contextualized to the benchmark of pooled-donor suspensions of PHHs (50-donor pool). Moreover, the extents of metabolic clearance (i.e., loss of parent compound) and capacity for formation of Phase II glucuronidation and sulfation metabolites were evaluated for a subset of Tox21 cell lines following in vitro exposures to 7-hydroxycoumarin (7HC), BPAF, and BPS. Given the purported lack of xenobiotic metabolism in cell lines used in Tox21 assays, our findings reveal both anticipated and surprising results with implications to the interpretation of Tox21 and in vitro toxicology bioactivity assay data. The study provides a benchmark to contextualize the proficiency of metabolic competence with in vitro model systems,

and highlights important considerations for chemicals with direct Phase II sulfation and glucuronidation metabolism.

2. Material and methods

2.1. Chemicals and reagents

All metabolism substrates and assay reagents were of the highest grade/purity (e.g., HPLC grade). Phenacetin, bupropion, midazolam (MI), acetaminophen, hydroxybupropion, 1-hydroxy midazolam, aminobenzotriazole (ABT), omeprazole (OMP), rifampicin (RIF), 7HC, 7-hydroxycoumarin glucuronide (7HCG), 7-hydroxycoumarin sulfate (7HCS), BPAF and BPS were obtained from Sigma-Aldrich (Saint Louis, MO, USA). Acetonitrile, methanol and formic acid were purchased from Fisher Scientific (Waltham, MA, USA). Eagle’s MEM medium (catalog No. 30–2003) was purchased from ATCC (Manassas, VA, USA) and was used for HepG2 and Hela cells culture. DMEM medium (catalog No. 11995) was purchased from Invitrogen (Carlsbad, CA, USA) and was used for HEK293 cells culture. McCoy’s 5a medium (catalog No. 30–2007) was purchased from ATCC (Manassas, VA, USA) and was used for HCT116 and ME-180 cells culture. F-12 K medium (catalog No. 30–2004) was purchased from ATCC (Manassas, VA, USA) and was used for CHO-K1 cells culture. DMEM/F12, GlutaMax™ supplement medium (catalog No. 10565) was purchased from Life Technologies (Carlsbad, CA, USA) and was used for GH3.TRE-Luc cells culture. Eagle’s Basal medium (catalog No. 21010) was purchased from Invitrogen (Carlsbad, CA, USA) and was used for C3H10T1/2 cells culture. MEM/EBSS medium (catalog No. SH 30024.01) was purchased from Hyclone (Seattle, WA, USA) and was used for MCF7 cells culture.

2.2. Cell lines, cell incubations and cell culture

2.2.1. Primary human hepatocytes

Cryopreserved pooled PHHs (from 50 individual donor preparations) were purchased from Life Technologies/Thermo (HUE50-D) and were stored in liquid nitrogen vapor until use. Hepatocytes were thawed using hepatocyte thawing medium (Lonza Biosciences Solutions; Morrisville, NC, USA). Cell viability and yield were determined using a Cellometer (Nexcelom Bioscience, Lawrence, MA, USA) with acridine orange/propidium iodide (AOPI) staining per the manufacturer’s protocol. Once thawed, hepatocyte maintenance medium (Lonza Biosciences Solutions; Morrisville, NC, USA) was added to attain a concentration of 1×10^6 cells/mL. 50 μ L of cell suspension was added to a V-bottom 96-well plate (Corning, NY, USA) at 50,000 cells/well with 50 μ L maintenance medium (control) or with 50 μ L of substrate at $2 \times$ final desired concentration in maintenance medium. The plates were incubated at 37 °C, 5% CO₂ with individual metabolism substrates (specific activity assays) or parent chemicals for metabolic clearance and metabolite formation assays (7HC, BPAF, BPS). PHHs were incubated for no more than 240 min in suspension form due to their loss of viability after 4 h.

2.2.2. Cell lines used for Tox21 chemicals screening (Tox21 cell lines)

Tox21 cell lines including HepG2, HEK293, Hela, HCT116, ME180, CHO-K1, GH3.TRE-Luc, C3H10T1/2 and MCF7 cells were obtained from National Center for Advancing Translational Sciences (NCATS, Rockville, MD) and cultured analogous to established cell culture conditions for each respective line. Evaluation of enzymatic specific activities, metabolic clearance, and metabolite formation was performed in Costar flat bottom 96-well cell culture plates (Corning, NY, USA) analogous to culture conditions (<https://tripod.nih.gov/tox21/assays/>) used for each respective Tox21 cell line (e.g., culture media and initial seeding densities). Prior to initiating metabolism incubations, cultures were allowed to reach approximately 80–90% confluence (Hsieh et al., 2017).

2.3. Cytochrome P450 specific activity assays

2.3.1. Metabolism incubations

PHHs and Tox21 cell lines (i.e., HepG2, HEK293, Hela, HCT116, ME180, CHO-K1, GH3.TRE-Luc, C3H10T1/2 and MCF7) were each incubated at 37 °C, 5% CO₂ with phenacetin, bupropion and midazolam (i.e., individually) with shaking (200 rpm) on an orbital shaker (model BT1502, Benchmark Scientific, Edison, NJ, USA). Phenacetin o-deethylase, bupropion hydroxylase, and midazolam 1-hydroxylation activities were assessed using *in situ* metabolism evaluations of metabolite formation in cell culture medium as previously described metabolic competence evaluations (Jackson et al., 2016; Ramaiahgari et al., 2017). All samples were stored at -80 °C prior to analysis by Liquid Chromatography-Mass Spectrometry (LC-MS).

2.3.2. CYP450 specific activities

The formation of acetaminophen, hydroxybupropion and 1-hydroxymidazolam was quantified by LC-MS. The analytical method used was a modification of protocol described earlier, and a complete list of instrumentation and analysis parameters are described in Supplementary Table 1–3. (Ramaiahgari et al., 2017). Briefly, 50 µL of incubation medium samples from each well were transferred to a 96-well plate and 170 µL of 0.1% formic acid in acetonitrile was added along with 10 µL of isotopically labeled internal standard mixture (Supplementary Table 1, 250 ng/mL each in acetonitrile). Samples were mixed, centrifuged (4000 rpm for 5 min) and 100 µL of the supernatant was transferred to a Waters Acquity 96-well plate and mixed with 100 µL of methanol:water (1:1). Samples were analyzed using an API 5000 triple quadrupole mass spectrometer (SCIEX, Concord, Ontario, Canada) coupled to a Waters Acquity HPLC system (Waters, Milford, MA, USA). A Phenomenex Luna C18 (50 × 2 mm, 5 µm) Column (Phenomenex, Torrance, CA, USA) with LC gradient method consisting of mobile phase A (0.05% formic acid in 5 mM ammonium formate) and mobile phase B (0.05% formic acid in acetonitrile:methanol (95:5)) at 2% mobile phase B increased linearly to 5% over 0.5 min, followed by Mobile phase B increase to 71% over 2 min and then return to 2% over 0.5 min. The injection volume was 10 µL. LC flow rates were 0.7 mL/min for acetaminophen method and 0.5 mL/min for hydroxybupropion and hydroxymidazolam method.

Standard curves ranged from 1 to 1000 ng/mL. Standards were prepared by mixing 50 µL of blank media with 10 µL of spiking solution for each compound at each respective concentration level and extracted as per standards. Accuracy of all standards were within 15% of nominal for all concentrations except at limit of quantitation (LOQ, 1 ng/mL). A quadratic equation with 1/x weighting was used to relate peak area response ratio of analyte to internal standard and concentration. Concentrations of specific metabolites were calculated by using the analyte peak area/internal standard peak area and using regression analysis in Analyst 1.6.2 (SCIEX, Concord, Ontario, Canada) to determine the ng/mL quantities. Metabolite formation rates were normalized to number of cells and incubation time, and expressed as pmol/min-million cells. Michaelis-Menten curve fitting was performed and graphed with GraphPad Prism (GraphPad Software, San Diego, CA).

2.4. Glucuronidation and sulfation activity

2.4.1. Probe substrates exposures

7HC (50 µM), BPAF (50 µM), and BPS (50 µM) were used as probe substrates to evaluate glucuronidation and sulfation proficiency with 3 cell lines used in Tox21 assays: HepG2, HEK293, and MCF7. Assay plates were incubated at 37 °C, 5% CO₂ with these drug substrates individually with shaking (200 rpm) on an orbital shaker (model BT1502, Benchmark Scientific, Edison, NJ, USA). Reactions were stopped at 150 min (PHH & Tox21 cell lines) and 24 h (Tox21 cell lines) by addition of 100 µL cold acetonitrile to each well, and plates were then centrifuged at 100 xg for 10 min. 100 µL aliquots of solvent-crashed lysates were subsequently transferred to analytical plates and stored in -80 °C until

analysis.

2.4.2. Glucuronidation and sulfation LC-MS

Phase II metabolites for 7HC, 7HCS and 7HCG were analyzed using a targeted Selected Ion Monitoring (SIM) method, while metabolites of BPAF and BPS were analyzed using a semi-targeted method with data collected in a full scan MS experiment but using inclusion lists of mass-to-charge values corresponding to glucuronide and sulfate conjugates of BPAF and BPS (Supplementary Table 4–5) using a Q Exactive Plus (Thermo Fisher, Waltham, MA, USA) with a Vanquish Ultra performance Liquid Chromatography (UPLC) system (Thermo Fisher, Waltham, MA, USA). Both methods used a Hypersil Gold aq Column (100 × 2.1 mm, 1.9 µm) (Thermo Fisher, Waltham, MA, USA). LC conditions consisted of a gradient method where mobile phase A was 0.1% formic acid in water and mobile phase B was 0.1% formic acid in acetonitrile. The gradient method started with a 2-min equilibration at 10% mobile phase B. Mobile phase B was then increased to 100% over 3 min and then held at 100% for 2 min. In the last minute of the gradient, mobile phase B was decreased back to 10%. LC flow rate was 0.3 mL/min. The injection volume was 2 µL. The column temperature was held at ambient temperature whereas the autosampler tray temperature was set at 4 °C. The mass spectrometer was calibrated with Thermo Fisher Scientific positive mode calibration solution. A complete list of MS analysis parameters is described in Supplementary Table 6–8.

2.5. Metabolic clearance assays

2.5.1. Chemicals exposure

Metabolic clearance assays following disappearance of parent compound over time were performed to determine respective rates of molecular transformation for 7HC, BSAF, and BPS in 50-donor pooled suspensions of PHH and 3 Tox21 cell lines (HepG2, HEK293, and MCF7) that produced detectable metabolite concentrations during specific activity evaluations. For these exposures, 1 µM 7HC, BPAF, or BPS prepared in cell culture medium. Assay plates were incubated at 37 °C, 5% CO₂ with individual substrates with shaking (200 rpm) on an orbital shaker (model BT1502, Benchmark Scientific, Edison, NJ, USA). Reactions were stopped at times 0, 15, 30, 60, 90, 120, 240, and 1440 min by adding 100 µL of cold acetonitrile to each well (i.e., whole well crash), and centrifuged at 100 xg for 10 min. 100 µL of each acetonitrile-crashed lysate were transferred to analytical plates and stored in -80 °C for subsequent LC-MS analysis.

2.5.2. Metabolic clearance LC-MS

Parent compounds 7HC, BPAF, and BPS were detected using a single ion monitoring method with MS/MS used for analyte confirmation. All data were collected using a Q Exactive Plus (Thermo Fisher, Waltham, MA, USA) with a Vanquish UPLC system (Thermo Fisher, Waltham, MA, USA). All compounds were analyzed in negative mode following LC separation with a Hypersil Gold aq Column (100 × 2.1 mm, 1.9 µm) (Thermo Fisher, Waltham, MA, USA) using a gradient method. LC conditions consisted of a gradient method where mobile phase A was 0.1% formic acid in water and mobile phase B was 0.1% formic acid in acetonitrile. The gradient method started with a 2-min equilibration at 10% mobile phase B. Mobile phase B was then increased to 30% over 2 min and then to 100% over the next 3 min. Gradient was held at 100% mobile phase B for a minute before being decreased back to 10% over the last minute of the run. The mass spectrometer was calibrated with Thermo Fisher Scientific positive mode calibration solution. A complete list of MS analysis parameters is described in Supplementary Table 9. Data was processed using Xcalibur 4.1 (Thermo Fisher, Waltham, MA, USA). Parent substrate concentration versus time plots were created using GraphPad Prism and fitted using non-linear regression analysis to determine half-life.

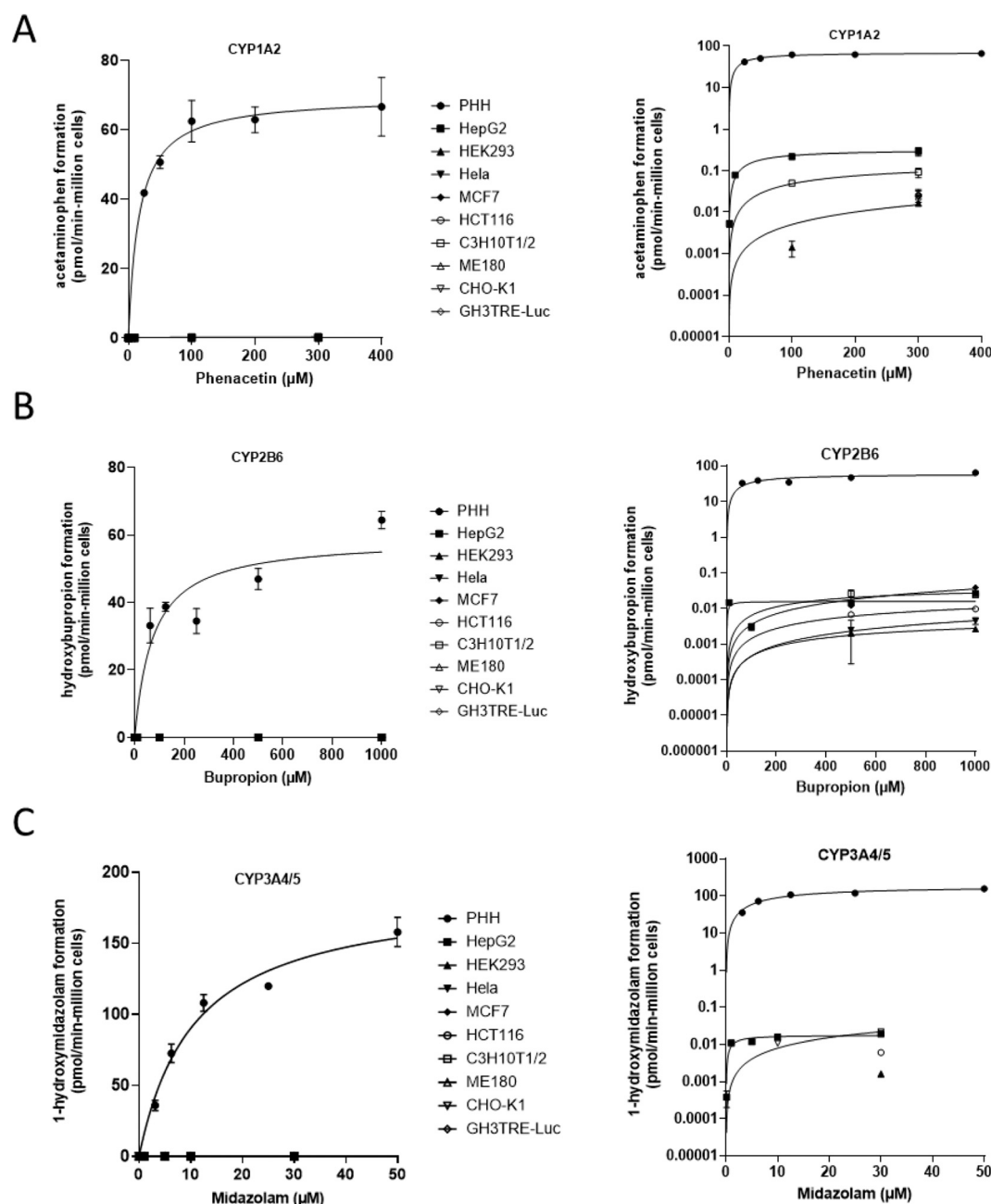


Fig. 1. Metabolism comparison between PHHs and Tox21 cell lines.

Fifty-donor pool of PHHs and Tox21 cell lines at $\sim 50,000$ cells/well were treated with probe substrates phenacetin (PHH = 25, 50, 100, 200, 400, Tox21 = 0.1, 1, 10, 100, 300) bupropion (PHH = 62.5, 125, 250, 500, 1000, Tox21 = 1, 10, 100, 500, 1000), and midazolam (PHH = 3.125, 6.25, 12.5, 25, 50, Tox21 = 0.1, 1, 5, 10, 30). LC-MS assays were performed to quantify acetaminophen (A: CYP1A2), hydroxybupropion (B: CYP2B6), and 1-hydroxymidazolam (C: CYP3A4/5/7) plotted on both linear (left) and log (right) scales. Metabolite concentrations were converted to pmol of total metabolite per million cells plated over the respective incubation periods used for each model (pmol/min- 10^6 cells).

2.6. CYP450-Glo bioluminescent assays

Luminescent screening assays were established with CYP450-Glo assays (Promega, Madison WI) for CYP1A2 and CYP3A4/5. These assays monitor the conversion by CYPs of inactive D-luciferin derivatives to an active form that produces light when firefly luciferase is added to the reaction mixture, with light intensity proportionally changing with enzymatic activity (Larson et al., 2011). Tox21 cell lines were cultured as a monolayer in 96-well white plate (Corning, NY, USA) at $\sim 50,000$ cells/well for 24 h. Then, cells were pretreated with the vehicle 0.1% DMSO (Basal), or OMP (100 μM), a well known CYP1A2 inducer, or ABT

(1 mM), a CYP1A2 inhibitor, for 24 h. In other experiments, cells were pretreated with the vehicle 0.1% DMSO (Basal), or RIF (20 μM), well known CYP3A4 inducer, or MI (30 μM), a CYP3A4 inhibitor, for 24 h. Assay background (BGCYP1A2 or BGCYP3A4) was measured in medium without cells. After treatment, the medium was replaced with fresh medium containing Luciferin-1A2 P450-Glo substrate (6 μM), or Luciferin-3A4 P450-Glo substrate (3 μM), cells were incubated for an additional 60 min. CYP450-Glo bioluminescent assays were performed by manufacturer's protocol using a CLARIOSTAR luminometer (BMG LABTECH).

Table 1Metabolism comparison between PHHs and Tox21 cell lines for acetaminophen formation (pmol/min-10⁶ cells); CYP1A2.

Phenacetin (μM)										
Cell Lines	Incubation Time	0.1	1	10	25	50	100	200	300	400
PHH	30 min	–	–	–	41.8 ± 0.353	50.7 ± 1.85	62.4 ± 6.00	62.8 ± 3.69	–	66.6 ± 8.43
HepG2	24 h	0.00	0.00527 ± 0.000925	0.0786 ± 0.00898	–	–	0.217 ± 0.0257	–	0.293 ± 0.0648	–
HEK293	24 h	0.00	0.00	0.00	–	–	0.00141 ± 0.000591	–	0.0163 ± 0.00150	–
Hela	24 h	0.00	0.00	0.00	–	–	0.00	–	0.00	–
MCF7	24 h	0.00	0.00	0.00	–	–	0.00	–	0.00	–
HCT116	24 h	0.00	0.00	0.00	–	–	0.00	–	0.0248 ± 0.00770	–
C3H10T1/2	24 h	0.00	0.00	0.00	–	–	0.0500 ± 0.00414	–	0.0914 ± 0.0236	–
ME180	24 h	0.00	0.00	0.00	–	–	0.00	–	0.0294 ± 0.00589	–
CHO-K1	24 h	0.00	0.00	0.00	–	–	0.00	–	0.0233 ± 0.00102	–
GH3TRE-Luc	24 h	0.00	0.00	0.00	–	–	0.00	–	0.0204 ± 0.00175	–

Table 2Metabolism Comparison between PHHs and Tox21 Cell Lines for Hydroxybupropion Formation (pmol/min-10⁶ cells); CYP2B6.

Bupropion (μM)										
Cell Lines	Incubation Time	0.1	1	10	62.5	100	125	250	500	1000
PHH	30 min	–	–	–	33.1 ± 5.14	–	38.7 ± 1.34	34.5 ± 3.70	46.9 ± 3.10	64.4 ± 2.60
HepG2	24 h	0.00	0.00	0.0146 ± 0.000928	–	0.00305 ± 0.000668	–	–	0.0141 ± 0.000381	0.0262 ± 0.00206
HEK293	24 h	–	0.00	0.00	–	0.00	–	–	0.00203 ± 0.000380	0.00262 ± 0.000174
Hela	24 h	–	0.00	0.00	–	0.00	–	–	0.00245 ± 0.000218	0.00447 ± 0.000947
MCF7	24 h	–	0.00	0.00	–	0.00	–	–	0.0119 ± 0.0206	0.0382 ± 0.00250
HCT116	24 h	–	0.00	0.00	–	0.00	–	–	0.00667 ± 0.000451	0.00958 ± 0.000485
C3H10T1/2	24 h	–	0.00	0.00	–	0.00	–	–	0.0257 ± 0.00337	0.0241 ± 0.000575
ME180	24 h	–	0.00	0.00	–	0.00	–	–	0.00	0.00
CHO-K1	24 h	–	0.00	0.00	–	0.00	–	–	0.00	0.00
GH3TRE-Luc	24 h	–	0.00	0.00	–	0.00	–	–	0.00	0.00

Table 3Metabolism Comparison between PHHs and Tox21 Cell Lines for 1-Hydroxymidazolam Formation (pmol/min-10⁶ cells); CYP3A4/5.

Midazolam (μM)										
Cell Lines	Incubation Time	0.1	1	3	5	10	12.5	25	30	50
PHH	30 min	–	–	35.8 ± 3.65	72.4 ± 6.53	–	107.9 ± 5.87	119.80 ± 1.41	–	157.9 ± 10.3
HepG2	24 h	0.00038 ± 0.00018	0.0109 ± 0.000694	–	0.0118 ± 0.000588	0.0157 ± 0.000372	–	–	0.0192 ± 0.000124	–
HEK293	24 h	–	0.00	–	0.00	0.00	–	–	0.00167 ± 0.000283	–
Hela	24 h	–	0.00	–	0.00	0.00	–	–	0.00	–
MCF7	24 h	–	0.00	–	0.00	0.00	–	–	0.00	–
HCT116	24 h	–	0.00	–	0.00	0.00	–	–	0.00607 ± 0.000435	–
C3H10T1/2	24 h	–	0.00	–	0.00	0.00	–	–	0.00	–
ME180	24 h	–	0.00	–	0.00	0.00	–	–	0.00	–
CHO-K1	24 h	–	0.00	–	0.00	0.0108 ± 0.00180	–	–	0.0217 ± 0.00212	–
GH3TRE-Luc	24 h	–	0.00	–	0.00	0.00	–	–	0.00	–

Table 4Specific activity maximum enzymatic velocity (V_{\max}) and Michelis-Menten equilibrium constant (K_m) estimates summary.

Cell Type	Incubation Time	V_{\max}			K_m		
		CYP1A2	CYP2B6	CYP3A4/5	CYP1A2	CYP2B6	CYP3A4/5
PHH	30 min	69.5 ± 2.12	59.0 ± 4.49	187 ± 8.35	16.7 ± 2.79	73.2 ± 23.3	10.7 ± 1.34
HepG2	24 h	0.320 ± 0.0234	0.0157 ± 0.00259	0.0177 ± 0.000912	38.4 ± 11.8	4.79 ± 5.48	0.969 ± 0.268
HEK293	24 h	NA	0.00616 ± 0.00221	NA	NA	1260 ± 734	NA
Hela	24 h	NA	NA	NA	NA	NA	NA
MCF7	24 h	NA	NA	NA	NA	NA	NA
HCT116	24 h	NA	0.0286 ± 0.0110	NA	NA	1910 ± 1050	NA
C3H10T1/2	24 h	0.168 ± 0.0393	0.0422 ± 0.0119	NA	247 ± 111	566 ± 354	NA
ME180	24 h	NA	NA	NA	NA	NA	NA
CHO-K1	24 h	NA	NA	0.175 ± 0.300	NA	NA	209 ± 405
GH3TRE-Luc	24 h	NA	NA	NA	NA	NA	NA

NA - Not applicable due to insufficient data for model fit.

2.7. Statistical analysis


All data represent the mean ± standard deviation (SD) of three or four independent experiments ($n = 3$ or 4), unless otherwise stated, using Microsoft Excel. Graphs were plotted using GraphPad Prism 7 (GraphPad software, San Diego, CA). Student's *t*-test or ANOVA with subsequent Dunnett's test were used as appropriate, with a $p < 0.05$

considered significant.

3. Results

3.1. Phase I metabolism evaluation

To characterize the perceived limitations for Tox21 cell lines to

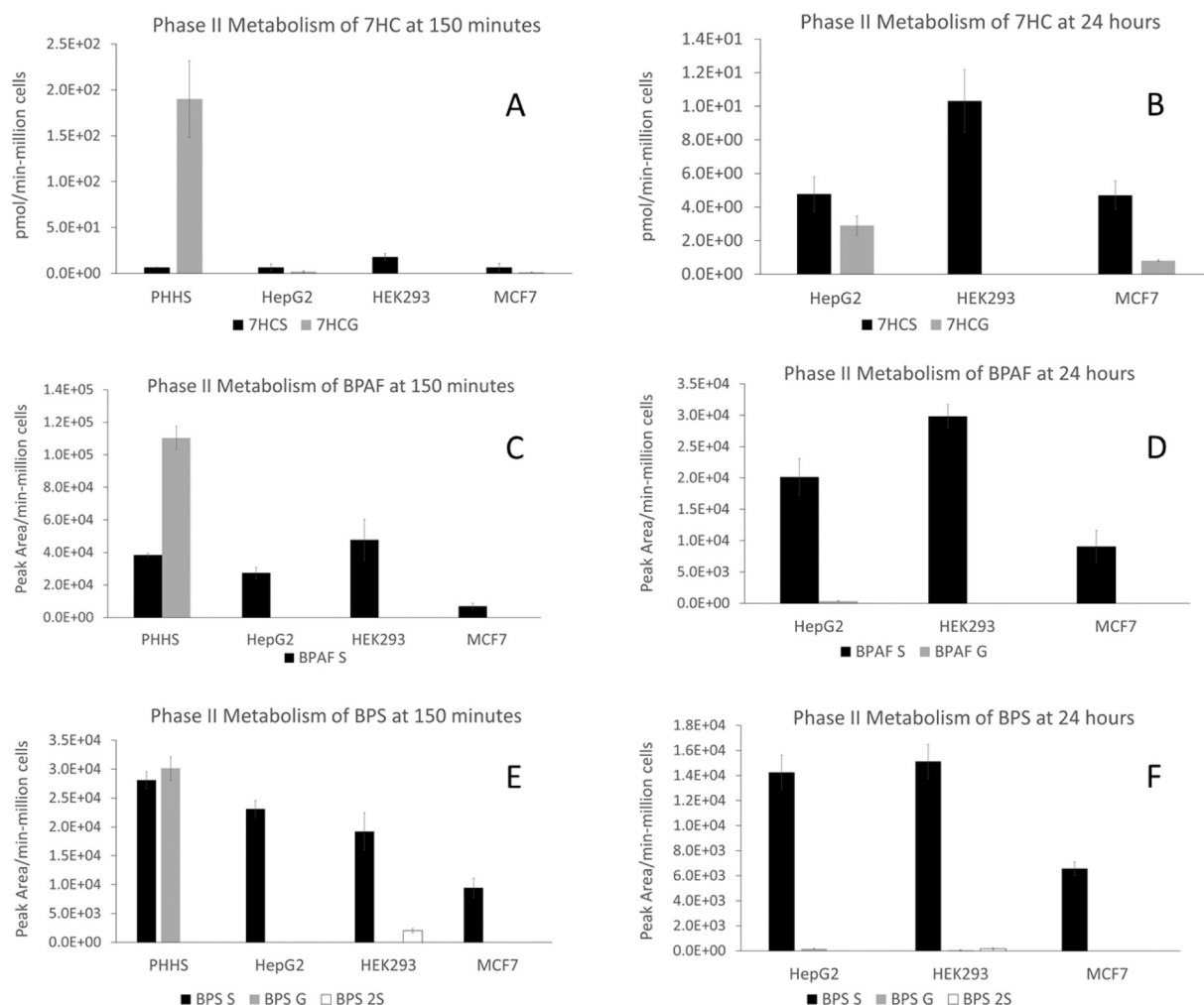**Fig. 2.** Glucuronide and sulfate conjugates assay in PHH and three Tox21 cell lines.

HepG2, HEK293 and MCF7 cells were incubated at 50,000 cells/well for 24 h individually in comparison to PHHs (50,000 cells/well) incubated with 50 μ M 7HC (A, 150 min; B, 24 h); 50 μ M BPAF (C, 150 min; D, 24 h), or 50 μ M BPS (E, 150 min; F, 24 h). LC-MS assays were performed to quantify 7HC-glucuronide and sulfate conjugates, semi-quantify BPAF-glucuronide, sulfate conjugates, BPS-glucuronide and sulfate conjugates. Metabolite concentrations were converted to pmol of total metabolite per million cells plated for 7HC conjugates (pmol/min- 10^6 cells) and peak area of total metabolite per million cells plated (peak area/min- 10^6 cells) for BPAF and BPS conjugates.

metabolize xenobiotics, HepG2, HEK293, Hela, HCT116, ME180, CHO-K1, GH3.TRE-Luc, C3H10T1/2 and MCF7 cultures were established and assayed in 96-well microplate formats for CYP specific activities to major human drug metabolizing enzymes CYP3A4/5, CYP1A2, and CYP2B6. Both substrate concentration and time course ranges were evaluated. As shown in Fig. 1 and Table 1-3, robust metabolite formation for marker metabolites with clinical drug substrates phenacetin, bupropion, and midazolam was observed after 30 min incubations with suspensions of PHHs, which appeared to plateau with increasing substrate concentrations consistent with Michaelis-Menten kinetics. In contrast, Tox21 cell lines cultured over a much longer 24 h incubation period produced profoundly lower specific activities requiring much longer incubation times to observe detectable metabolite formation. Determinable V_{max} and K_m values are shown in Table 4.

In general, hepatic lineage HepG2 cultures were the most proficient Tox21 cell line when summing CYP1A, CYP2B and CYP3A specific activities. The ranking for all Tox21 cell lines was HepG2 > C3H10T1/2 > CHO-K1 > MCF7 ~ HCT116 ~ ME180 ~ HEK293 ~ GH3TRE-Luc > Hela. For CYP1A metabolism, seven cell lines produced detectable acetaminophen metabolite levels following phenacetin substrate exposures (24 h). The relative ranking of activity was: HepG2 > C3H10T1/2 > ME180 > HCT116 > CHO-K1 > GH3TRE-Luc > HEK293 (Table 1).

For CYP2B activity, assessed by bupropion hydroxylation formation, 7 cell lines produced some detectable, but near the lower limit of quantitation activity after very long (i.e., 24 h incubation periods) to accumulate detectable metabolite. The ranking of CYP2B activities observed was MCF7 > HepG2 > C3H10T1/2 > HCT116 > Hela > HEK293 at 1000 μ M bupropion concentrations (Table 2). For CYP3A activity, only four cell line cultures produced detectable 1-hydroxymidazolam after 24 h incubations with the following ranking CHO-K1 ~ HepG2 > HCT116 > HEK293 (Table 3). All Tox21 culture model activities, when benchmarked to PHH suspensions, produced very low CYP3A-type metabolism activities that are widely considered the most important enzymes for drug/xenobiotic metabolism in humans. Hepatic lineage HepG2 were only weakly effective in metabolizing phenacetin to acetaminophen (<1% of PHHs), bupropion to hydroxybupropion (<0.1% of PHHs) and midazolam to 1-hydroxymidazolam (< 0.1% of PHHs). These data benchmark and define the highly limited capacity of Tox21 cell lines to model major human xenobiotic metabolism functionality, and indicate extremely limited capacity for CYP1A-, CYP2B-, or CYP3A-dependent metabolic transformation or bioactivation in response to chemical exposures.

To further evaluate the extent of xenobiotic metabolism functionality with Tox21 cell lines, inhibitors of xenobiotic metabolism were applied. As shown in Supplemental Fig. 1, HepG2, MCF7, CHO-K1, C3H10T1/2, and GH3TRE-Luc were assayed using CYP1A2 luminescent substrate. Comparison of normalized luminescence data showed MCF7 > HepG2 > CHO-K1 > C3H10T1/2 ~ GH3TRE-Luc, and marginally inhibitable metabolite formation after exposure to the P450 inhibitor ABT only inhibited in CHO-K1 cells. These data were consistent with the very low enzymatic activities observed, and further confirm their deficiency for CYP450 metabolism.

While it is generally known that immortalized cell lines lack proficient signal transduction pathways related to drug metabolizing enzymes (e.g., PXR, CAR), we further evaluated Tox21 cell line exposures with prototype inducers of xenobiotic metabolism. AhR activator OMP only marginal increased luminescence ranking: MCF7 > CHO-K1 > HepG2 > C3H10T1/2 > GH3TRE-Luc. Overall, CYP1A2/AhR functionality was highest in MCF7 among the Tox21 assays evaluated. For CYP3A/PXR functionality, HepG2 were superior in basal CYP3A luminescence with a ranking of HepG2 > MCF7 > HCT116 > C3H10T1/2 ~ Hela (Supplemental Fig. 2), which was appreciably inhibitable in all 5 cell lines with substrate-depleting high concentrations of MI with a ranking HepG2 > C3H10T1/2 > HCT116 > MCF7 > Hela. As

Table 5
Phase II metabolite formation rates.

150 min												
Metabolites	7HCS	7HCG	BPAF S	BPAF G	BPS S	BPS G	BPS 2S					
Cell Lines	pmol/min-million cells		Peak area/min-million cells		Peak area/min-million cells							
PHHs	6.55 ± 0.365	190 ± 41.8	38,500 ± 1110	110,000 ± 7240	28,100 ± 1440	30,100 ± 2010	1					
HepG2	6.46 ± 3.63	1.96 ± 0.887	27,400 ± 3450	1	23,100 ± 1470	1	1					
HEK293	17.9 ± 3.56	1	47,800 ± 12,600	1	19,200 ± 3210	1	1990 ± 378					
MCF7	6.63 ± 4.56	1.34 ± 0.206	7000 ± 1880	1	9460 ± 1680	1	1					
24 h												
Metabolites	7HCS	7HCG	BPAF S	BPAF G	BPS S	BPS G	BPS 2S					
Cell Lines	pmol/min-million cells		Peak area/min-million cells		Peak area/min-million cells							
PHHs	4.77 ± 1.04	2.89 ± 0.563	20,200 ± 2930	375 ± 48.0	14,300 ± 1380	170 ± 11.4	2					
HepG2	10.3 ± 1.86	1	29,800 ± 1830	1	15,100 ± 1360	53.8 ± 6.89	1					
HEK293	4.69 ± 0.848	0.809 ± 0.0457	9090 ± 2570	1	6560 ± 527	1990 ± 378	1					
MCF7												

¹ Metabolite not detected.

² PHH suspension cells no longer viable at this time point.

Table 6
Chemicals clearance comparison between PHHs and Tox21 cell lines.

Chemicals	Concentration (μ M)	Cell Culture Model	Percent Depletion 240 min (% \pm SD)	Percent Depletion 24 h (% \pm SD)	Mean Half-life (min)*	95% Confidence Interval Half-life (min)	Fit Method Description
7-HC	1	PHHs (50-donor pool)	99 \pm 5.9	Not determined	22.7	19.8–26.1	One-phase exponential decay
7-HC	1	HepG2	27 \pm 6.8	76 \pm 11	290.6	196–512	One-phase exponential decay
7-HC	1	HEK293	50 \pm 12	[†] 99 \pm 12	209.9	173–259	One-phase exponential decay
7-HC	1	MCF7	40 \pm 18	99 \pm 12	292.5	196–530	One-phase exponential decay
BPAF	1	PHHs (50-donor pool)	82 \pm 2.9	Not determined	86.3	51.2–202	One-phase exponential decay
BPAF	1	HepG2	36 \pm 8.2	66 \pm 6.6	31.2 (fast) 1626 (slow)	very wide	Two-phase exponential decay
BPAF	1	HEK293	36 \pm 5.4	60 \pm 4.5	43.0 (fast) ~1887 (slow)	very wide	Two-phase exponential decay
BPAF	1	MCF7	42 \pm 7.7	55 \pm 7.4	23.3 (fast) ~9684 (slow)	very wide	Two-phase exponential decay
BPS	1	PHHs (50-donor pool)	59 \pm 6.7	Not determined	205.2	73.5–268	One-phase exponential decay
BPS	1	HepG2	17 \pm 11	39 \pm 7.4	326.1	164–3053	One-phase exponential decay
BPS	1	HEK293	9.5 \pm 11	36 \pm 12	NC	NC	One-phase exponential decay
BPS	1	MCF7	22 \pm 7.0	26 \pm 7.3	119.7	48.9–272	One-phase exponential decay

NC: Not calculated due to insufficient trend within dataset for fit model estimates.

* value from 150 min.

[†] value applies lower limit of quantification (LLOQ) where mean concentrations less than LLOQ.

expected, no response to prototype PXR agonist rifampicin (RIF) was observed in any of the Tox21 cell lines evaluated, with the exception of a small increase in CYP3A4 activity in Hela cells.

3.2. Phase II metabolism evaluation

Phase II metabolism in Tox21 cell lines was evaluated with bisphenols BPAF and BPS, along with 7-HC, and metabolite formations were assayed using LC-MS. Direct overlay of existing Tox21 assay data for estrogen receptor (ER) activation revealed BPAF to be more potent than BPA or BPS (Supplemental Fig. 3), thus investigation of their in vitro metabolism was explored (<https://sandbox.ntp.niehs.nih.gov/tox21-curve-visualization> (Huang et al., 2014)). Preliminary experiments revealed detectable concentrations of 7HC-glucuronide (7HCG) and sulfate (7HCS) only in HepG2, HEK293 and MCF7 (data not shown). Based on these data, more extensive evaluation of glucuronide and sulfate metabolites of 7HC, BPAF, BPS (50 μ M for each) was evaluated over 150 min and 24 h incubations (i.e., 7HCG, 7HCS, BPAFG, BPAFC, BPSG, and BPSS). As shown in Fig. 2 and Table 5, PHH suspensions produced far more 7HCG than 7HCS conjugate (Fig. 2A) than Tox21 cell lines after 150 min. Since PHH suspensions lose cell viability after 2–4 h in culture, it was not possible to evaluate the 24 h time point. After 24 h, HepG2 and MCF7 cultures achieved substantial amounts of 7HCG, whereas HEK293 produced no detectable 7HCG at either 150 min or 24 h. However, 7HCS conjugates were detected in all 3 Tox21 cell lines (i.e., HepG2, HEK293 and MCF7) after both 150 min and 24 h, demonstrating a greater proficiency for sulfation over glucuronidation.

BPAFG conjugate was readily observed in PHH suspensions after 150-min incubations, with higher proportions of BPAFS vs. BPAFG (Fig. 2C). By comparison, very low BPAFG was observed in these 3 Tox21 cell lines. Interestingly, BPAF sulfation was clearly observed in all 3 Tox21 cell lines after 150-min incubations, with levels rivaling those produced in PHH suspensions. Only HepG2 cultures retained observable levels of both sulfate- and glucuronide-conjugated BPAF after 24 h (Fig. 2D), consistent with its hepatic lineage.

As shown in Fig. 2E & 2F, BPS incubations produced robust BPSS and

BPSG metabolites equivalent to 150-min PHH incubations. In contrast, no detectable BPS G was observed in Tox21 cell lines after 150-min incubations. BPSS was observed at comparable levels to PHH suspensions, consistent with results with 7HC and BPAF. After 24 h incubations, BPSS concentrations further increased for all 3 Tox21 cell lines and was analogous to the profile observed for BPAFS with proportionally small amounts of BPSG. BPS2S conjugates (i.e., 2 sulfate conjugates) were also assayed and detectable in HEK293 cells after both 150 min and 24 h (Fig. 2E, F). Taken together, these data demonstrate that Tox21 cell lines HepG2, HEK293, and MCF7 do display robust Phase II sulfation capacity with all three substrates evaluated (7HC, BPAF, & BPS), yet minimally supported Phase II glucuronidation.

3.3. Metabolic clearance evaluation

In vivo, the spectrum of adaptive and toxicological responses to xenobiotic exposure occurs as a result of the combined dynamic exposure to both parent chemical and downstream metabolites over time. Therefore, to better understand the impact of observed metabolite formations, the extent of metabolic clearance (i.e., loss of parent form of chemical to alternate forms) was further evaluated with HepG2, HEK293, and MCF7 cell lines. For this, 7HC, BPAF, and BPS were assayed over 24 h (Table 6, Figs. 3–5). Analogous to Phase II metabolite formation evaluations, 7HC was measured via LC-MS to monitor loss of 7HC (1 μ M initial substrate concentration) and derive percent depletion and apparent half-life ($T_{1/2}$) values (Table 6). 7HC was extensively ‘cleared’ (i.e., 99% conversion to other chemical structure forms) in PHH suspensions after 240-min incubations, which revealed an estimated half-life of ~23 min (Table 6, Fig. 3). Twenty-four-hour incubations with PHH suspensions were not possible due to their short-term viability (i.e., ~2–4 h). In contrast, HepG2 cells revealed 27% depletion after 240 min with an estimated half-life of ~290 min. Interestingly, HEK293 had a faster loss of 7HC relative to hepatic lineage HepG2 with 50% depletion after 240 min, and 99% turnover after 24 h with an overall $T_{1/2}$ ~210 min. This observation was consistent with the elevated extent of 7HC sulfation observed in Fig. 2. MCF7 produced a

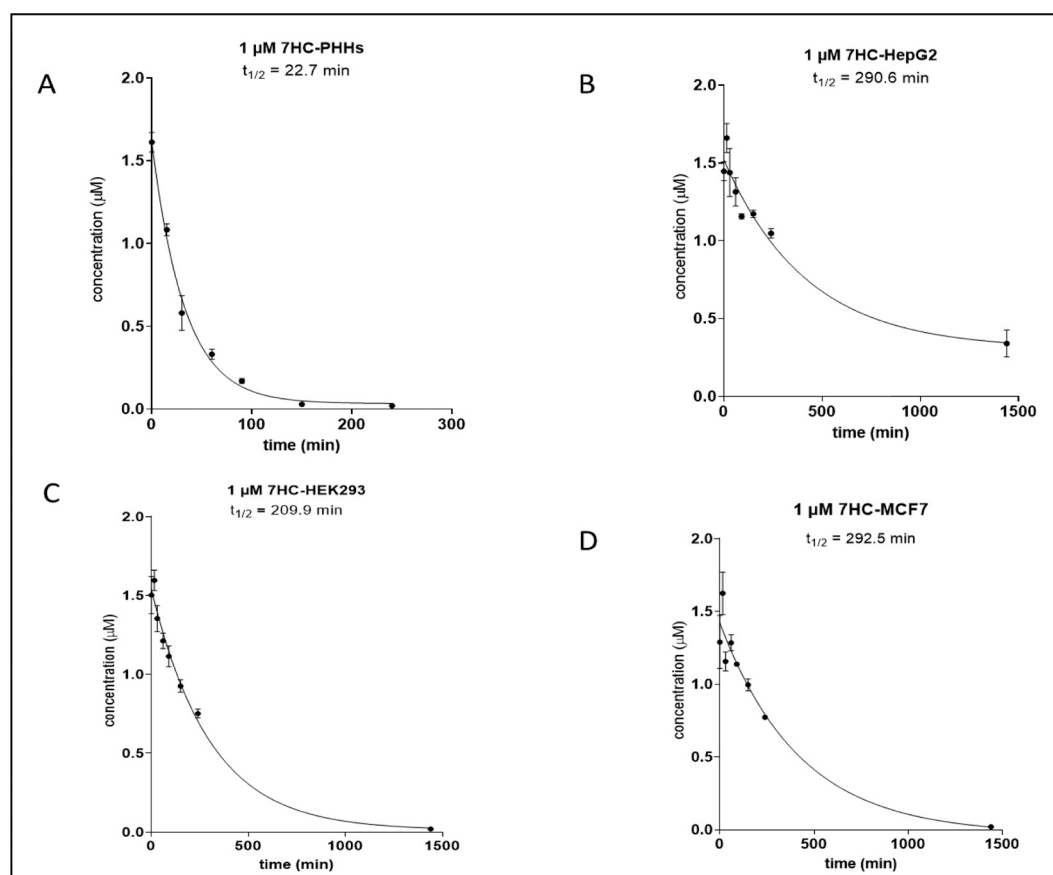


Fig. 3. Clearance of 7HC by PHH and three Tox21 cell lines.

HepG2, HEK293 and MCF7 cells were incubated at 50,000 cells/well for 24 h individually in comparison to PHHs (50,000 cells/well) incubated with 1.0 μM 7HC at 0, 15, 30, 60, 90, 150 and 240 min. LC-MS analysis was performed to monitor the clearance of 7HC parent compound at 1.0 μM to derive percent depletion and apparent half-life ($t_{1/2}$). 7HC clearance in PHHs (A); 7HC clearance in HepG2 cells (B); 7HC clearance in HEK293 cells (C); 7HC clearance in MCF7 cells (D).

comparable extent of 7HC depletion with 40% depletion after 240 min with $T_{1/2}$ –290 min (Fig. 3D).

The metabolic clearance of BPAF was extensive with 7HC in PHH suspensions with 82% depletion and $T_{1/2}$ –86 min (Table 6, Fig. 4). HepG2 cells revealed an apparent biphasic depletion with initial fast and delayed slow disappearance of BPAF with $T_{1/2}$ –31 and 1600 min, respectively. Percent depletion of BPAF in HepG2 cells was 36% for 240 min and 66% after 24 h. This extensive depletion of parent BPAF demonstrates Tox21 cell lines are in fact competent for xenobiotic metabolism by conversion parent chemicals to alternative but form, likely to be sulfation based on our observations of Phase II metabolite data. This non-zero parent chemical depletion has the potential to impact the interpretation of Tox21 assay data. HEK293 showed a similar biphasic depletion of BPAF with $T_{1/2}$ –43 and 1900 min, respectively. MCF7 further mirrored the biphasic depletion of BPAF revealing the fastest initial metabolic clearance rate of across these three Tox21 cell lines with $T_{1/2}$ –23 min, and the slowest second phase with $T_{1/2}$ –9700 min. These data were consistent with the extensive BPAF sulfate formation shown in Fig. 2.

Metabolic clearance was by far the least extensive for BPS among these 3 compounds (i.e., 7HC, BPAF, and BPS). As shown in Table 6 and Fig. 5, depletion of BPS was 59% after 240 min in PHH suspensions with correspondingly slower $T_{1/2}$ –205 min. In HepG2, significantly lower depletion of 17% was observed with BPS after the initial 240 min, and 39% turnover of BPS after 24 h with overall half-life $T_{1/2}$ –326 min. Compared with HepG2, BPS turnover in HEK293 was only 9.5% after 240 min and 36% after 24 h and an overall half-life that was not calculated due to insufficient trend within the dataset. By comparison, MCF7 produced the most extensive BPS depletion among Tox21 cell

lines with ~30% depletion, which was largely achieved in the initial 240 min and ranked as the fastest BPS depletion half-life of $T_{1/2}$ –120 min. These observations were largely consistent with metabolite formation data.

4. Discussion

The Tox21 high throughput screening program has been criticized for a perceived lack of ‘metabolic competence’ across the suite of bioactivity assays screening thousands of environmental and pharmaceutical chemicals. However, to date this speculative critique has not been experimentally confirmed or characterized. Given the importance of biotransformation and internal exposure profiles to toxicological outcomes, the perceived yet poorly characterized lack of xenobiotic metabolism in Tox21 cell culture models, and the broad impact of Tox21 assay data made available to the public for thousands of chemicals, it is important to characterize and contextualize the capacity of Tox21 cell culture models for ‘metabolic competence’. In this study, we established a framework to benchmark ‘metabolic competence’ using the hepatic ‘gold standard’ of pooled-donor PHH suspensions, and characterized 9 commonly used Tox21 cell culture models for major drug metabolizing enzyme activities (e.g., CYP3A, CYP1A, CYP2B, UGT, SULT). The results highlight both anticipated and surprising results, which have implications to the field of in vitro toxicology and in vitro to in vivo translation of bioactivity assay data.

As we considered how best to approach the challenges in characterizing ‘metabolic competence’ within cell culture models used for Tox21 assays, we focused on CYP, UGT, and SULT metabolism assays. The rationale focused on their preeminent role in xenobiotic metabolism

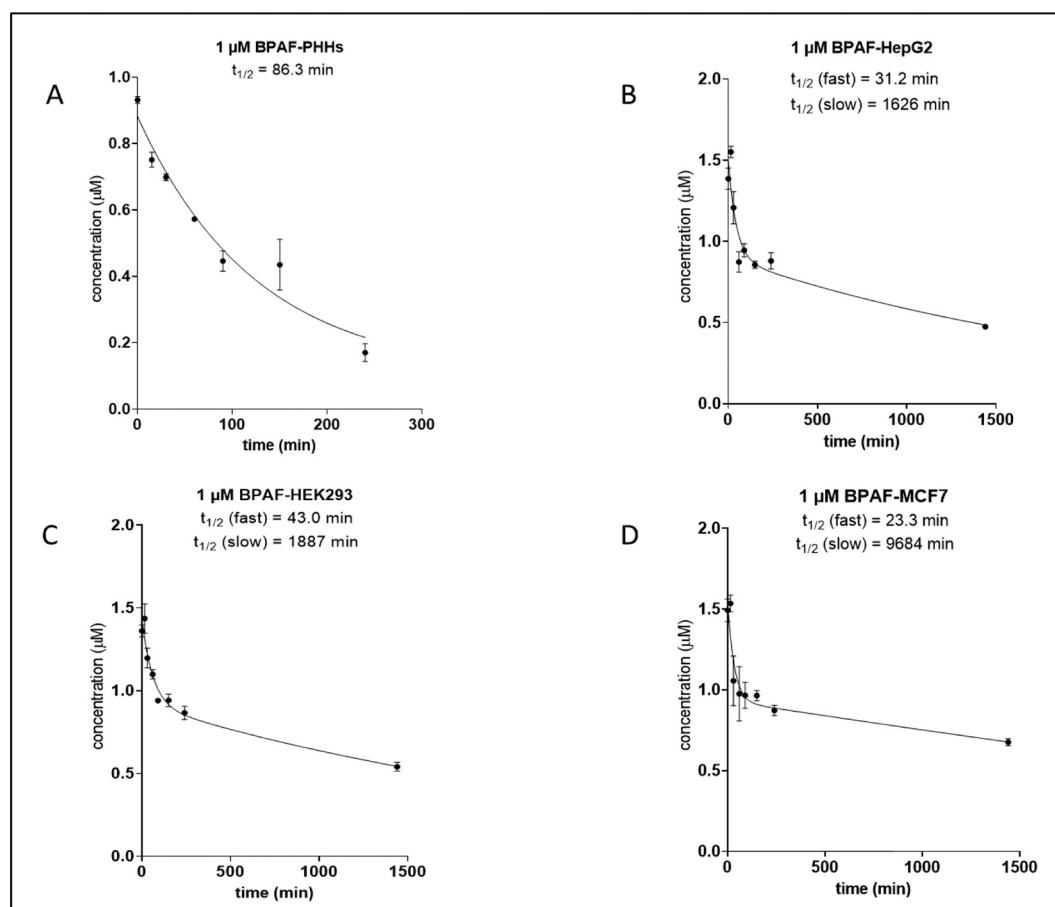


Fig. 4. Clearance of BPAF by PHH and three Tox21 cell lines.

HepG2, HEK293 and MCF7 cells were incubated at 50,000 cells/well for 24 h individually in comparison to PHHs (50,000 cells/well) incubated with 1.0 μ M BPAF at 0, 15, 30, 60, 90, 150 and 240 min. LC-MS analysis was performed to monitor the clearance of BPAF parent compound at 1.0 μ M to derive apparent half-life ($t_{1/2}$). BPAF clearance in PHHs (A); BPAF clearance in HepG2 cells (B); BPAF clearance in HEK293 cells (C); BPAF clearance in MCF7 cells (D).

and clearance in humans, and our need to address the predominant modes of metabolic activation of chemicals to toxic metabolites that historically involve hepatic metabolism (Walsh, 2006). Xenobiotic metabolism is known to both detoxify and activate toxicological responses, and can arise from both hepatic and extrahepatic metabolism (Jhajra et al., 2012). The enzymes commonly involved in metabolic activation often include UGTs (e.g., acyl glucuronide electrophiles), peroxidases, monoamine oxidases, and most often cytochromes P450 (CYPs). CYPs form a variety of reactive metabolites including epoxides, quinones, heterocyclics, halocarbons, hydrazines, and free radical ions, and are attributed to bioactivation roles with well-known toxicants such as benzo(a)pyrene, aflatoxin B1, benzene, dimethylnitrosamine, carbon tetrachloride, and acetaminophen. In humans, CYPs drive the majority of xenobiotic metabolism clearance via CYP3A, CYP1A, CYP2C, enzymes, and to a less prevalent extent via CYP2B, CYP2E, and CYP2D enzymes. These enzymes typically initiate Phase I metabolism to transform lipophilic chemicals towards more hydrophilic forms poised for Phase II metabolism and excretion. This rationale guided our approach to evaluate, major CYP, UGT, and SULT metabolite formations, and metabolic clearance in context with PHHs suspensions. PHH suspensions integrate active/passive uptake transport, intrinsic production of requisite enzyme cofactors at physiological proportions (e.g., NADPH, UDPGA, PAPS), and functional expression of a physiologically-relevant spectrum of xenobiotic metabolizing enzymes produced in human liver within the initial 1–2 h of incubation. The differentiation state of PHHs in culture can be inferred based on their CYP450 activities that dynamically change based on culture conditions. We observed that

Tox21 cell lines were largely devoid of reference CYP450 metabolite formation relative to PHH suspensions using established clinical drug substrates for CYP3A4/5, CYP1A2, and CYP2B6 as sentinels to broader ‘metabolic competence’ (Jackson et al., 2016; Lu et al., 2015). It is important to note that while these specific activities (i.e., metabolite formation of acetaminophen, hydroxybupropion, and 1-hydroxymidazolam) are well-established for human liver enzymes, their specificity for metabolizing enzymes with non-human origin Tox21 cell lines CHO-K1 (hamster), GH3TRE (rat), and C3H3T1/2 (mouse) are likely more generalized reflections of CYP3A, CYP1A, and CYP2B metabolism as some conservation of metabolism is known across a CYP subfamily. Another consideration to contextualize these findings is that many of these Tox21 cell lines are not of hepatic lineage, thus lower activities for major xenobiotic metabolizing enzymes were anticipated given the importance of hepatic metabolism in drug clearance. Extrahepatic human tissues often vary in the degree and proportions of CYP expression, which can impact the specificity of metabolite formation due to substrate overlap among enzymes. A final point of context is that CYPs involved in endogenous metabolism such as CYP19A1 (i.e., aromatase catalyzing conversion of androgens to estrogens) with high specificity for endogenous substrates and narrow ranges of substrate recognition are most likely operative in Tox21 culture models, but generally do not serve as rate-determining xenobiotic metabolism/clearance. These considerations guided our rationale to evaluate a broad spectrum of xenobiotic metabolism substrate concentrations over shorter and longer (i.e., 24 h) incubations with Tox21 cultures. Overall, these data definitively establish a profound lack of CYP-dependent ‘xenobiotic

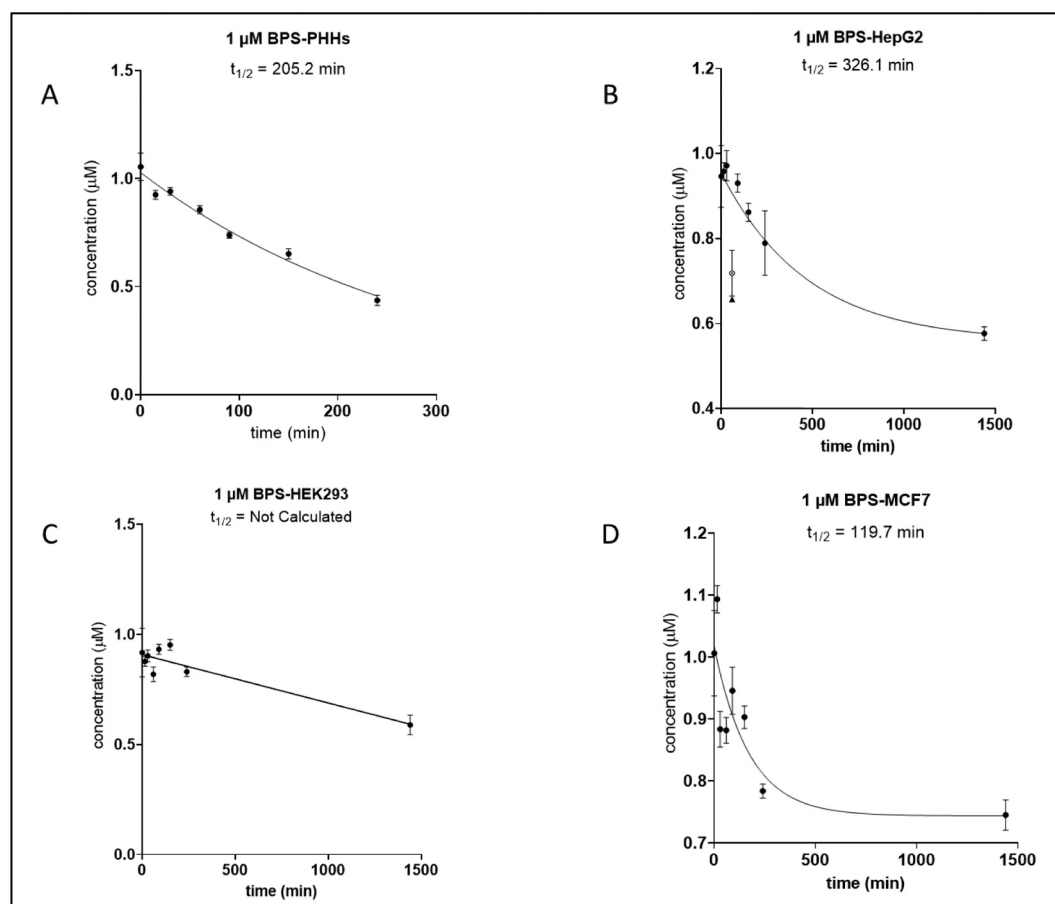


Fig. 5. Clearance of BPS by PHH and three Tox21 cell lines.

HepG2, HEK293 and MCF7 cells were incubated at 50,000 cells/well for 24 h individually in comparison to PHHs (50,000 cells/well) incubated with 1.0 μM BPS at 0, 15, 30, 60, 90, 150 and 240 min. LC-MS analysis was performed to monitor the clearance of BPS parent compound at 1.0 μM to derive apparent half-life ($t_{1/2}$). BPS clearance in PHHs (A); BPS clearance in HepG2 cells (B); BPS clearance in HEK293 cells (C); BPS clearance in MCF7 cells (D).

metabolism', as defined within the rationale of this study, and establish a benchmarking framework to evaluate CYP-related 'metabolic competence' with in vitro toxicology models.

The most surprising finding within this study was the substantial capacity of Tox21 cell culture models to produce Phase II sulfate metabolites of 7HC, BPAF, and BPS that in some instances rivaled PHH suspensions. These data were further supported by metabolic clearance data indicating significant loss of parent chemical over 240 min and 24 h exposures. These data demonstrated a clear 'competence' to appreciably biotransform phenolic chemicals. These data also showed a strong bias towards sulfation in Tox21 cell lines compared to PHH suspensions with more balanced glucuronidation and sulfation. The extents of metabolic clearance (i.e., disappearance of parent chemical over time) mirrored the extents of Phase II metabolite formation to sulfation/glucuronidation in most data sets. These data suggest that Phase II metabolism is driving the loss of parent BPS, BPAF, and 7HC in these cultures, but further study with inhibitors towards a mass balance evaluation would be required to confirm these coincident observations. It is interesting to note that with Tox21 cultures, though proficient for sulfation of BPAF and BPS akin to PHH suspensions, appeared to reach a plateau of parent chemical depletion with 1 μM exposures across multiple cell lines. This could reflect additional rate-determining factors such as free fraction protein binding (e.g., 10% FBS) or active transport processes that compartmentalize parent chemicals and limit the extent of turnover. In addition, UGT and SULT enzymes often display substrate affinity K_m values in the higher micromolar to low millimolar range, which may also play a role in the observed limit of further BPS and BPAF depletion between 240 min and 24 h. For context, Ishii and Tamura reported the

apparent K_m of bisphenol A for sulfotransferase activity was 74 μM , and the depletion of bisphenol A was only 4% after 5 days incubations with 10 μM exposure (Ishii and Tamura, 2003). Given the range of exposures and shorter (< 24 h in most assays) exposure times for Tox21 assay systems, we could hypothesize from these observations that some phenolic and aromatic hydroxyl compounds, along with other sulfotransferase substrates, may deplete to some extent during Tox21 assays. Future analysis is required to assess the potential impact of these findings across the thousands of Tox21 chemicals evaluated to date. Depletion, along with the potential for skewed metabolite exposure profiles favoring sulfation, could result in shifted potencies that under- or over-estimate the potential for human biological activity and parent chemical responses in humans. These data highlight the challenge in ascribing Tox21 assay results to nominal exposure concentrations of parent chemicals that presume a lack of 'metabolic competence'. Furthermore, a more nuanced definition of 'metabolic competence' reflective of tissue-specific biotransformation is needed. As we consider the fit-for-purpose use of different types of in vitro liver models, both HepG2 and PHHs likely have a useful role to fulfill. For high throughput screening, HepG2 are a practical and inexpensive model system with some degree of CYP1A metabolism, and surprisingly robust sulfation capacity. However, their profoundly limited CYP450 metabolism and hepatic receptor signaling pose fundamental challenges for translational hepatotoxicity predictions. Suspension cultures of PHHs remain the 'gold standard' for accurate estimations of human metabolism outcomes, but their extremely limited viability time frame (i.e., a few hours due to the absence of epithelial-like cell-cell and cell matrix interactions) essentially reflects the responses of dying cells that are incapable of

toxicological modeling without more advanced culture configurations (e.g., sandwich cultures, 3D spheroids, microphysiological systems). Future studies to evaluate the in vitro disposition of Tox21 compounds should further improve and refine our understanding of the impacts of in vitro exposure profiles and 'skewed' metabolism on the translation of assay results to humans.

In conclusion, we have characterized and benchmarked a profound deficiency of major human CYP450 drug metabolism activities with 9 cell line backgrounds used in Tox21 assays relative to 50-donor pooled suspensions of primary human hepatocytes. Thus, compounds requiring hepatic CYP450 bioactivation would likely be poorly modeled with existing Tox21 assay systems. This approach establishes a framework for evaluation of human 'metabolic competence' with in vitro systems using pooled-donor suspensions of PHHs. Focused evaluation of phenolic and aromatic hydroxyl-containing compounds 7HC, BPAF, and BPS towards Phase II metabolism revealed a surprisingly robust capacity for sulfation over glucuronidation metabolism, which were consistent with metabolic clearance assay data. These results indicate a more nuanced definition of 'metabolic competence' is required as we expand the use of Tox21 assay data, and these findings likely have implications for the interpretation and quantitative translation of a subset of chemicals within existing Tox21 assay data appreciably metabolized by sulfotransferase enzymes.

Declaration of Competing Interest

Authors declare that they have no conflicts of interest that might be relevant to the content of the manuscript. The authors received no specific funding for this work.

Acknowledgments

We thank Drs. Esra Mutlu and Nisha Sipes for their helpful reviews of this manuscript. We want to thank Dr. Raymond Tice and John Bucher for their guidance in exploring these xenobiotic metabolism challenges for the Tox21 Program, and also thank Jui-Hua Hsieh for her assistance with Tox21 ER agonist data.

Quantitation of CYP activity assays was performed by RTI (research Triangle Park, NC) for the National Toxicology Program, National Institute of Environmental Health Sciences, National Institutes of Health, U.S. Department of Health and Human Services, under contract HHSN273201400022C.

NIH IRP project ES103318-04 (2019); Biomolecular Screening and Alternative Approaches for the National Toxicology Program; <https://intramural.nih.gov/search/index.taf>.

Appendix A. Supplementary data

Supplementary data to this article can be found online at <https://doi.org/10.1016/j.tiv.2021.105109>.

References

- Collins, F.S., Gray, G.M., Bucher, J.R., 2008. Toxicology. Transforming environmental health protection. *Science* 319, 906–907. <https://doi.org/10.1126/science.1154619>.
- Extrahepatic Drug-Metabolizing Enzymes and Their Significance. In: Jhagra, Shalu, Ramesh Varkhede, Ninad, Ahire, Deepak Suresh, Naik, Bukke Vidyasagar, Prasad, Bhagwat, Paliwal, Jyoti, Singh, Saranjit (Eds.), October 2012. Wiley, 10.1002/9780470921920.edm028.
- Fischer, F.C., Abele, C., Henneberger, L., Klüber, N., König, M., Mühlenbrink, M., Schlichting, R., Escher, B.I., 2020. Cellular metabolism in high-throughput in vitro reporter gene assays and implications for the quantitative in vitro-in vivo extrapolation. *Chem. Res. Toxicol.* 33, 1770–1779. <https://doi.org/10.1021/acs.chemrestox.0c00037>.
- Gerin, B., Dell'Aiera, S., Richert, L., Smith, S., Chanteux, H., 2013. Assessment of cytochrome P450 (1A2, 2B6, 2C9 and 3A4) induction in cryopreserved human hepatocytes cultured in 48-well plates using the cocktail strategy. *Xenobiotica* 43, 320–335. <https://doi.org/10.3109/00498254.2012.719088>.
- Guengerich, F.P., 1997. Role of cytochrome P450 enzymes in drug-drug interactions. *Adv. Pharmacol.* 43, 7–35. [https://doi.org/10.1016/s0009-2797\(97\)00068-9](https://doi.org/10.1016/s0009-2797(97)00068-9).
- Hewitt, N.J., Lecluyse, E.L., Ferguson, S.S., 2007. Induction of hepatic cytochrome P450 enzymes: methods, mechanisms, recommendations, and in vitro-in vivo correlations. *Xenobiotica* 37, 1196–1224. <https://doi.org/10.1080/00498250701534893>.
- Hsieh, J.H., Huang, R., Lin, J.A., Sedykh, A., Zhao, J., Tice, R.R., Paules, R.S., Xia, M., Auerbach, S.S., 2017. Correction: real-time cell toxicity profiling of Tox21 10K compounds reveals cytotoxicity dependent toxicity pathway linkage. *PLoS One* 12, e0181291. <https://doi.org/10.1371/journal.pone.0181291>.
- Huang, R., Sakamuru, S., Martin, M.T., Reif, D.M., Judson, R.S., Houck, K.A., Casey, W., Hsieh, J.H., Shockley, K.R., Ceger, P., Postel, J., Witt, K.L., Tong, W., Rotroff, D.M., Zhao, T., Shinn, P., Simeonov, A., Dix, D.J., Austin, C.P., Kavlock, R.J., Tice, R.R., Xia, M., 2014. Profiling of the Tox21 10K compound library for agonists and antagonists of the estrogen receptor alpha signaling pathway. *Sci. Rep.* 4, 5664–5672. <https://doi.org/10.1038/srep05664>.
- Ishii, Y., Tamura, H.O., 2003. Sulfoconjugation of bisphenol A in a human neuroblastoma cell line, NB-1. *J. Health Sci.* 49, 319–323. <https://doi.org/10.1248/jhs.49.319>.
- Jackson, J.P., Li, L., Chamberlain, E.D., Wang, H., Ferguson, S.S., 2016. Contextualizing hepatocyte functionality of cryopreserved HepaRG cell cultures. *Drug Metab. Dispos.* 44, 1463–1479. <https://doi.org/10.1124/dmd.116.069831>.
- Kafert-Kasting, S., Alexandrova, K., Barthold, M., Laube, B., Friedrich, G., Arseniev, L., Hengstler, J.G., 2006. Enzyme induction in cryopreserved human hepatocyte cultures. *Toxicology* 220, 117–125. <https://doi.org/10.1016/j.tox.2005.12.013>.
- Kavlock, R.J., Austin, C.P., Tice, R.R., 2009. Toxicity testing in the 21st century: implications for human health risk assessment. *Risk Anal.* 29, 485–497. <https://doi.org/10.1111/j.1539-6924.2008.01168.x>.
- Larson, B., Moeller, T., Banks, P., Cali, J.J., 2011. Automated triplexed hepatocyte-based viability and CYP1A and -3A induction assays. *J. Biomol. Screen.* 16, 895–902. <https://doi.org/10.1177/1087057111411482>.
- Lee, K.S., Kim, S.K., 2013. Direct and metabolism-dependent cytochrome P450 inhibition assays for evaluating drug-drug interactions. *J. Appl. Toxicol.* 33, 100–108. <https://doi.org/10.1002/jat.1720>.
- Liao, C., Liu, F., Guo, Y., Moon, H.B., Nakata, H., Wu, Q., Kannan, K., 2012. Occurrence of eight bisphenol analogues in indoor dust from the United States and several Asian countries: implications for human exposure. *Environ. Sci. Technol.* 46, 9138–9145. <https://doi.org/10.1021/es302004w>.
- Lu, J.T., Einhorn, S., Venkatarangan, L., Miller, M., Mann, D.A., Watkins, P.B., LeCluyse, E., 2015. Morphological and functional characterization and assessment of iPSC-derived hepatocytes for in vitro toxicity testing. *Toxicol. Sci.* 147, 39–54. <https://doi.org/10.1093/toxsci/kfv117>.
- Mao, J., Tay, S., Khojasteh, C.S., Chen, Y., Hop, C.E., Kenny, J.R., 2016. Evaluation of time dependent inhibition assays for marketed oncology drugs: comparison of human hepatocytes and liver Microsomes in the presence and absence of human plasma. *Pharm. Res.* 33, 1204–1219. <https://doi.org/10.1007/s11095-016-1865-9>.
- Ramaiahgari, S.C., Waidyanatha, S., Dixon, D., DeVito, M.J., Paules, R.S., Ferguson, S.S., 2017. From the cover: three-dimensional (3D) HepaRG spheroid model with physiologically relevant xenobiotic metabolism competence and hepatocyte functionality for liver toxicity screening. *Toxicol. Sci.* 159, 124–136. <https://doi.org/10.1093/toxsci/kfx122>.
- Song, S., Ruan, T., Wang, T., Liu, R., Jiang, G., 2012. Distribution and preliminary exposure assessment of bisphenol AF (BPAF) in various environmental matrices around a manufacturing plant in China. *Environ. Sci. Technol.* 46, 13136–13143. <https://doi.org/10.1021/es303960k>.
- Tice, R.R., Austin, C.P., Kavlock, R.J., Bucher, J.R., 2013. Improving the human hazard characterization of chemicals: a Tox21 update. *Environ. Health Perspect.* 121, 756–765. <https://doi.org/10.1289/ehp.1205784>.
- Walsh, J.S., 2006. Metabolic activation-role in toxicity and idiosyncratic reactions. In: Borchardt, R.T., Kerns, E.H., Hageman, M.J., Thakker, D.R., Stevens, J.L. (Eds.), *Optimizing the "Drug-Like" Properties of Leads in Drug Discovery. Biotechnology: Pharmaceutical Aspects, vol IV.* Springer, New York, NY.
- Wang, L., Zhang, Z., Xu, X., Zhang, D., Wang, F., Zhang, L., 2015. Simultaneous determination of four trace level endocrine disrupting compounds in environmental samples by solid-phase microextraction coupled with HPLC. *Talanta* 142, 97–103. <https://doi.org/10.1016/j.talanta.2015.04.043>.
- Wetmore, B.A., Wambaugh, J.F., Ferguson, S.S., Sochaski, M.A., Rotroff, D.M., Freeman, K., Clewell 3rd, H.J., Dix, D.J., Andersen, M.E., Houck, K.A., Allen, B., Judson, R.S., Singh, R., Kavlock, R.J., Richard, A.M., Thomas, R.S., 2012. Integration of dosimetry, exposure, and high-throughput screening data in chemical toxicity assessment. *Toxicol. Sci.* 125, 157–174. <https://doi.org/10.1093/toxsci/kfr254>.
- Yamasaki, K., Takeyoshi, M., Yakabe, Y., Sawaki, M., Takatsuki, M., 2003. Comparison of the reporter gene assay for ER-alpha antagonists with the immature rat uterotropic assay of 10 chemicals. *Toxicol. Lett.* 142, 119–131. [https://doi.org/10.1016/s0378-4274\(03\)00019-5](https://doi.org/10.1016/s0378-4274(03)00019-5).